


HOSPITAL CARE

Surflo™

Winged Infusion Set Product Information

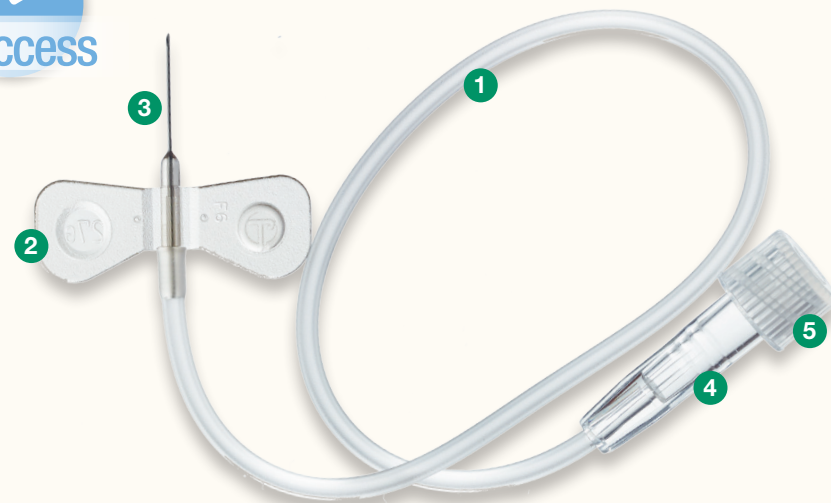

The wings of the Surflo™ Winged Infusion Sets are designed to ensure stability and grip during venipuncture. Thanks to the slightly angled wings (< 3° in relation to the cannula), less pressure is put on the blood vessel when puncturing.

Description	Raw Materials
1 Tube	Polyvinyl chloride
2 Wing hub	Polyvinyl chloride + coloured master batch
3 Cannula	Stainless steel
4 Adapter	PMMA
5 Adapter cap	Polypropylene
Protector	Polyethylene

Does not contain components made of DEHP or natural rubber latex.

Sterilization

With ethylene oxide to a Sterility Assurance Level of at least 10⁻⁶

Storage

- Shelf life of 5 years („NL“ items) or 3 years („NLSA“ & „ELF“ items)
- The lot number and expiry date are printed on the unit pack
- Do not store at extreme temperature and humidity

Ease of use

The comfort of your patients is in your hands

Product Range

Product Code	Colour	External Cannula Ø		Cannula Length	Dead Space Volume	Tubing Length	Replacement for
SV*18NL30	Pink	18G	1.20 mm	19 mm	0.45 ml	300 mm	SV*18BLK03
SV*19NL30	Ivory	19G	1.10 mm	19 mm	0.45 ml	300 mm	-
SV*21NL30	Green	21G	0.80 mm	19 mm	0.41 ml	300 mm	SV*21BLK03
SV*22NL30	Black	22G	0.70 mm	19 mm	0.41 ml	300 mm	SV*22BLK03
SV*23NL30	Blue	23G	0.60 mm	19 mm	0.40 ml	300 mm	-
SV*25NL30	Orange	25G	0.50 mm	19 mm	0.40 ml	300 mm	-
SV*25NLSA02	Orange	25G	0.50 mm	9.5 mm	0.19 ml	90 mm	SV*25BLSA02
SV27ELF	Grey	27G	0.40 mm	13 mm	0.30 ml	200 mm	SV*27EL



Packaging

Unit Box	Shipping Carton		
	Quantity	Gross Weight	Dimensions (L x W x H)
50	500	3.2 kg	510 x 257 x 260 mm
50	500	3.1 kg	512 x 262 x 264 mm
50	500	3.2 kg	510 x 257 x 260 mm
50	500	3.2 kg	510 x 257 x 260 mm
50	500	3.1 kg	512 x 262 x 264 mm
50	500	3.0 kg	512 x 262 x 264 mm
50	500	2.7 kg	510 x 257 x 260 mm
50	500	2.8 kg	510 x 257 x 260 mm

Single unit: Blister pack (thermoformed blister film and paper)

Regulatory Status

CE Certification – Class IIa medical device

Notified Body TÜV SÜD 0123
Manufacturer TERUMO MEDICAL PRODUCTS HANGZHOU,
 M4-9-5, Economic & Technological
 Development Zone, Hangzhou, 310018 China

EC Representative TERUMO EUROPE N.V.
 Interleuvenlaan 40, 3001 Leuven, Belgium

Safety warnings and (pre)cautions

- Do not use if package is damaged
- Do not reuse
- Dispose of safely after single use according to local regulations to avoid risk of infection.

Always consult the product label for a complete overview of warnings, cautions and/or precautions prior to actual use.

Choose from a wide range from 18G up to 27G

The outside diameter of the needle can be recognized by the ISO colour-coded wing hub.

Your Terumo partner



TERUMO EUROPE N.V.

Terumo Corporation **Terumo Europe N.V.**
 +81 3 6742 8500 +32 16 38 12 11

www.terumo-europe.com

All brand names are trademarks or registered trademarks owned by TERUMO CORPORATION, its affiliates, or unrelated third parties.
 Published by Terumo Europe N.V.

The image is indicative. Colours, shapes, sizes, components may differ in real life.