

Title: Labelling Specification for Exufiber Ag

Page 1(7)

POUCH

Making Sure that the Correct Top Web and reference are correct for the job / size of dressing being produced


REF	SKU	Range	Size	Top Web to be used
603400	603400-01	CE	2x45cm	417400-xx
603401	603401-01	CE	5x5cm	417401-xx
603402	603402-01	CE	10x10cm	417402-xx
603403	603403-01	CE	15x15cm	417403-xx
603404	603404-01	CE	4.5x10cm	417404-xx
603405	603405-01	CE	4.5x20cm	417405-xx
603406	603406-01	CE	4.5x30cm	417406-xx
603407	603407-01	CE	20x30cm	417407-xx
603412	603412-00	CE	10x10cm	417402-xx
603420	603420-00	US	0.8x17.7" (2x45cm)	417420-xx
603421	603421-00	US	2x2" (5x5cm)	417421-xx
603422	603422-00	US	4x4.8" (10x12cm)	417422-xx
603423	603423-00	US	6x6" (15x15cm)	417423-xx
603424	603424-00	US	8"x12"(20x30cm)	417424-xx
603425	603425-00	US	10x10cm	417425-xx

xx – Released version to be used

A 2d matrix containing the GTIN, Lot number, and expiry date (as per PD-517213) should be printed in line on the pouch, this information should also be applied in a human readable format along with the product reference. CE marked products also have the MD symbol and a Sterile barrier symbol printed in line.

See example of printing and placement below. Printing should be clearly legible and free from error, barcodes should register a grade C or above on the Axicon verifier.


CE Mark variant (Items sold into the UK)

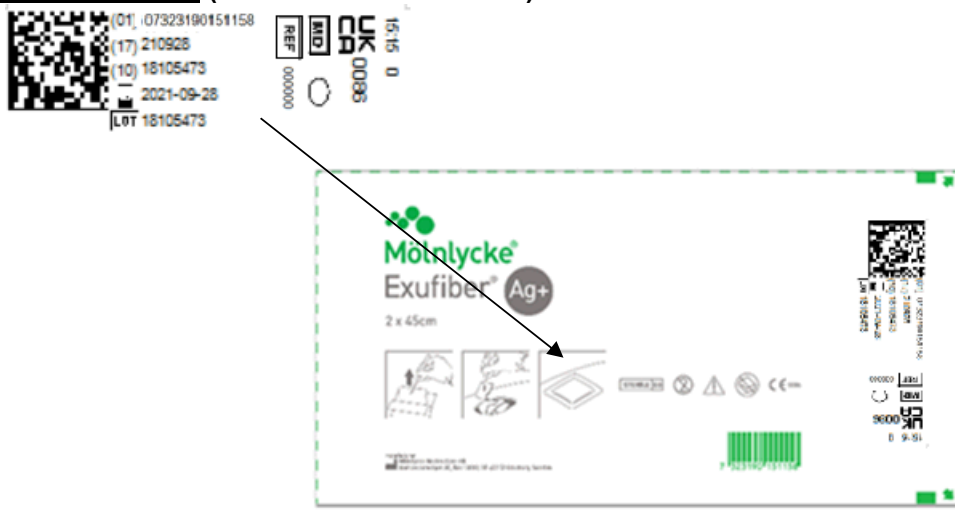


Title: Labelling Specification for Exufiber Ag

Page 2(7)

Spanish Only Variant

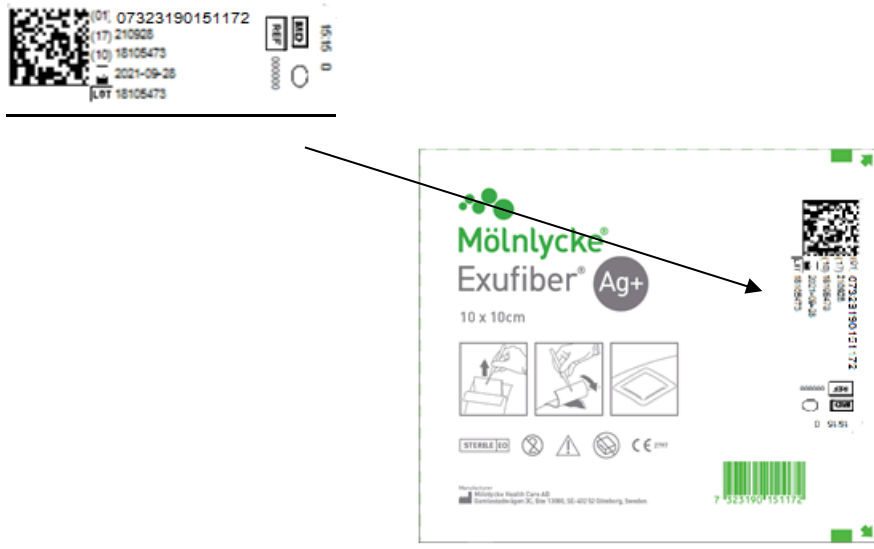

US Variant

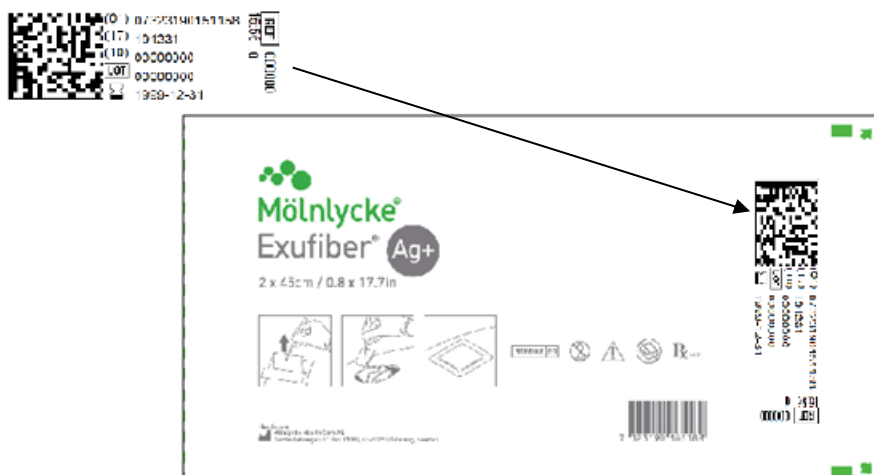

RETAILER

Making Sure that the Carton reference is correct for the job / size of dressing being produced

REF	SKU	Range	Size	Retail Carton	Reimbursement Label
603400	603400-01	CE	2x45cm	447400-xx	N/A
603401	603401-01	CE	5x5cm	447401-xx	N/A
603402	603402-01	CE	10x10cm	447402-xx	N/A
603403	603403-01	CE	15x15cm	447403-xx	N/A
603404	603404-01	CE	4.5x10cm	447404-xx	N/A
603405	603405-01	CE	4.5x20cm	447405-xx	N/A

Internal—Verify revision is current before using

Form ID: 000032, rev.04

Site(s): GLO

Title: Labelling Specification for Exufiber Ag


Page 3(7)

603406	603406-01	CE	4.5x30cm	447406-xx	N/A
603407	603407-01	CE	20x30cm	447407-xx	N/A
603412	603412-00	CE	10x10cm	447412-xx	496712-xx
603420	603420-00	US	0.8x17.7" (2x45cm)	447420-xx	N/A
603421	603421-00	US	2x2" (5x5cm)	447421-xx	N/A
603422	603422-00	US	4x4.8" (10x12cm)	447422-xx	N/A
603423	603423-00	US	6x6" (15x15cm)	447423-xx	N/A
603424	603424-00	US	8"x12"(20x30cm)	447424-xx	N/A
603425	603425-00	US	10x10cm	447425-xx	N/A

xx – Released version to be used

EU Variant

Example of 2d GTIN, Lot number and Expiry date and MD symbol placement in the varnish free area on the back of the retail box, printing should be clearly legible and free from error, barcodes should register a grade C or above on the Axicon verifier.

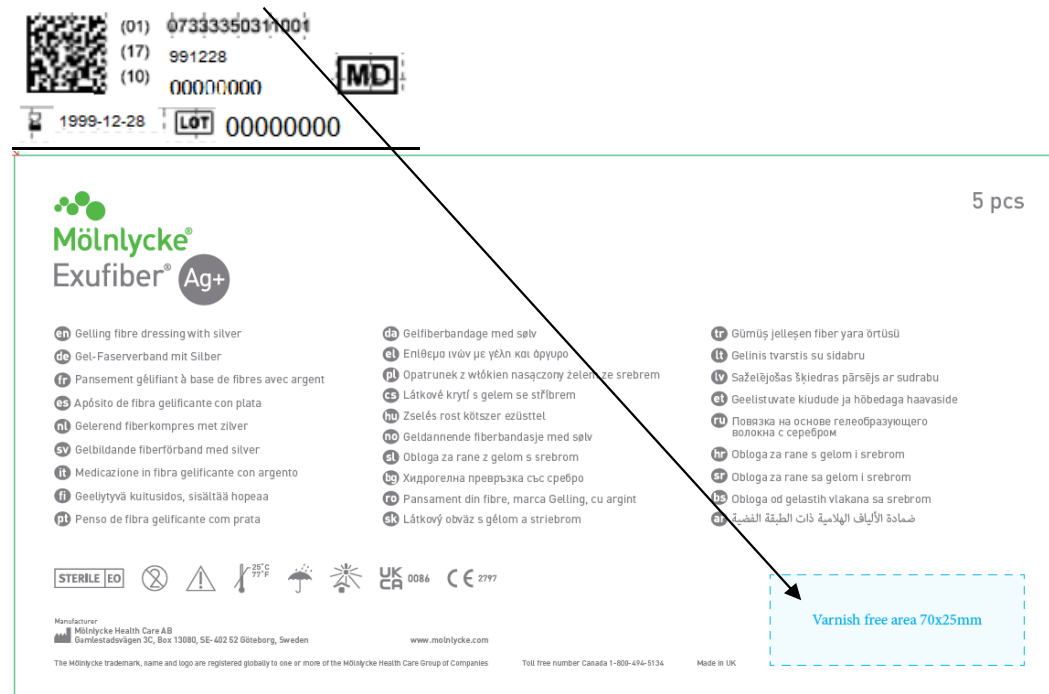


The diagram illustrates the back of the retail box for Mölnlycke Exufiber Ag+. It shows the following elements:

- Top Left:** A QR code and a 2D GTIN (0733335031001) with its corresponding numbers (01), (17), and (10).
- Top Center:** An MD symbol (MD) and a lot number (LOT) 00000000.
- Top Right:** An expiry date (1999-12-28).
- Bottom Left:** The Mölnlycke logo and product name "Exufiber Ag+".
- Bottom Center:** A list of multilingual descriptions of the product, including "Gelling fibre dressing with silver", "Gel-Faserverband mit Silber", "Pansement géifiant à base de fibres avec argent", "Apósito de fibra gelificante con plata", "Gelerend fiberkompres met zilver", "Gelbildande fiberförband med silver", "Medicazione in fibra gelificante con argento", "Geelytyvä kuitusidos, sisältää hopeaa", "Penso de fibra gelificante com prata", "Gel-fiberbandage med sølv", "Enlëþu ivón me yelín kai órpuro", "Opatrunek z włókien nasączony żelazem i srebrem", "Látkové krytí s gelem se stříbrem", "Zselés rost kötszer ezüsttel", "Geldannende fiberbandasje med sølv", "Obloga za rane z gelom i srebrom", "Хидрогелна превръзка със сребро", "Pansament din fibre, marca Gelling, cu argint", "Látkový obváz s gélom a striebrom", "Gümüş jelleyen fiber yara örtüsü", "Gelinis tvarstis su sidabru", "Saželajošas šķiedras pārsējs ar sudrabu", "Geelstuvate kiudude ja hõbedaga haavaside", "Повязка на основе гелеобразующего волокна с серебром", "Obloga za rane s gelom i srebrom", "Obloga za rane sa gelom i srebrom", "Obloga od gelastih vlakana sa srebrom", and "ضمادة الألياف الهلامية ذات الطبقة الفضية".
- Bottom Right:** A dashed blue box labeled "Varnish free area 70x25mm".
- Bottom Center:** A row of symbols including "STERILE EO", a crossed-out circle, a triangle, a thermometer, a sun, a star, and "UK 0086 CE 2797".
- Bottom Left (Small Text):** Manufacturer information: Mölnlycke Health Care AB, Gamlestadsvägen 3C, Box 13080, SE-402 52 Göteborg, Sweden. Website: www.molnlycke.com.
- Bottom Center (Small Text):** The Mölnlycke trademark, name and logo are registered globally to one or more of the Mölnlycke Health Care Group of Companies. Toll free number Canada 1-800-494-5134. Made in UK.

Title: Labelling Specification for Exufiber Ag


Page 4(7)

The Spanish Retail carton must also have a reimbursement label printed on the authorised blue label applied, placement illustrated below

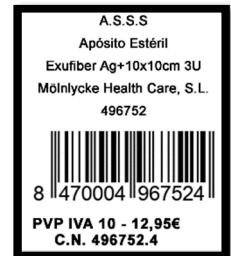

US Variant

Example of 2d GTIN, Lot number and Expiry date placement in the varnish free area on the back of the retail box, printing should be clearly legible and free from error, barcodes should register a grade C or above on the Axicon verifier.

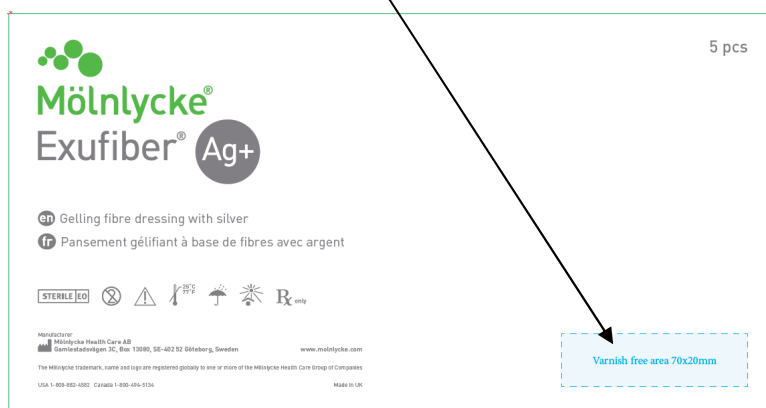

Title: Labelling Specification for Exufiber Ag

Page 5(7)

TRP

Making Sure that the kitting list is correct for the job / size of dressing being produced, scan the kitting list and the following Static information will be retrieved from our ERP system (Syspro), Reference, Re-order number, dressing size, Qty Breakdown, Number of pieces, and GS1 code the variable LOT and Expiry date is also populated from Syspro.

All sizes of dressing use the same blank label 307395-00



REF	SKU	Range	Dressing Size	Qty Breakdown	Number of Pieces per Shipper	GS1 Code
603400	603400-01	CE	2x45cm	(4 x 5)	20	07333350782221
603401	603401-01	CE	5x5cm	(4 x 10)	40	07333350981914
603402	603402-01	CE	10x10cm	(6 x 10)	60	07333350265557
603403	603403-01	CE	15x15cm	(6 x 10)	60	07333350129286
603404	603404-01	CE	4.5x10cm	(4 x 10)	40	07333350041830
603405	603405-01	CE	4.5x20cm	(5 x 10)	50	07333350763169
603406	603406-01	CE	4.5x30cm	(6 x 10)	60	07333350010492
603407	603407-01	CE	20x30cm	(4 x 5)	20	07333350240677
603412	603412-00	CE	10x10cm	(6 x 3)	18	07333350354701
603420	603420-00	US	0.8x17.7" (2x45cm)	(4 x 5)	20	07323190152476
603421	603421-00	US	2x2" (5x5cm)	(4 x 10)	40	07323190152551
603422	603422-00	US	4x4.8" (10x12cm)	(6 x 10)	60	07323190152599
603423	603423-00	US	6x6" (15x15cm)	(6 x 10)	60	07323190152674
603424	603424-00	US	8"x12"(20x30cm)	(4 x 5)	20	07323190152711
603425	603425-00	US	10x10cm	(6 x 10)	60	07323190152759

Title: Labelling Specification for Exufiber Ag


Page 6(7)

The TRP label should be attached to the short side of the box as illustrated. Printing should be clearly legible and free from error with the barcode achieving a grade C or above. Particular attention to check Reference, lot and expiry information are accurate.

Example of CE TRP Label Highlighted yellow information retrieved from Syspro.

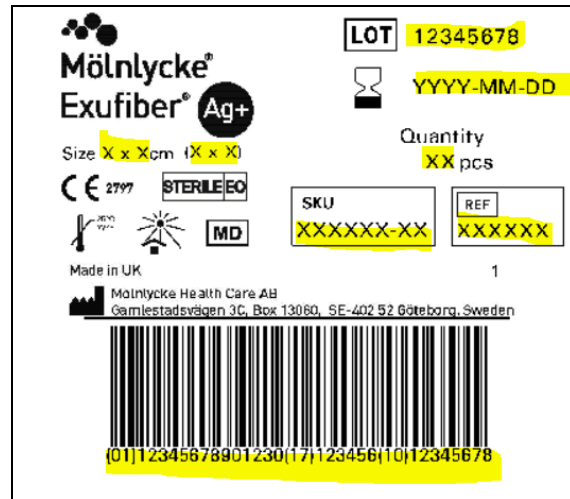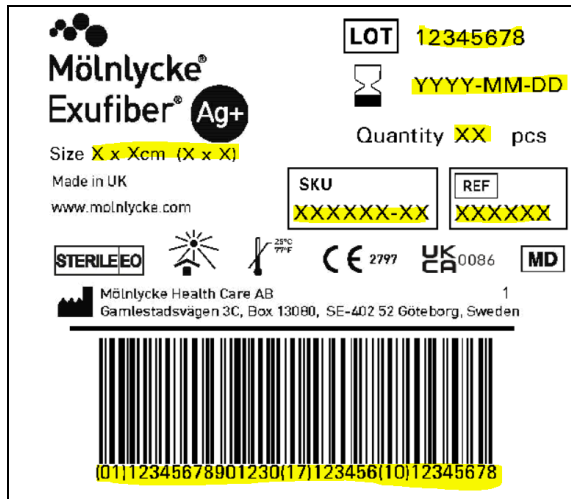Items sold into the UK
Spanish Only Variant


Example of US TRP Label Highlighted yellow information retrieved from Syspro.

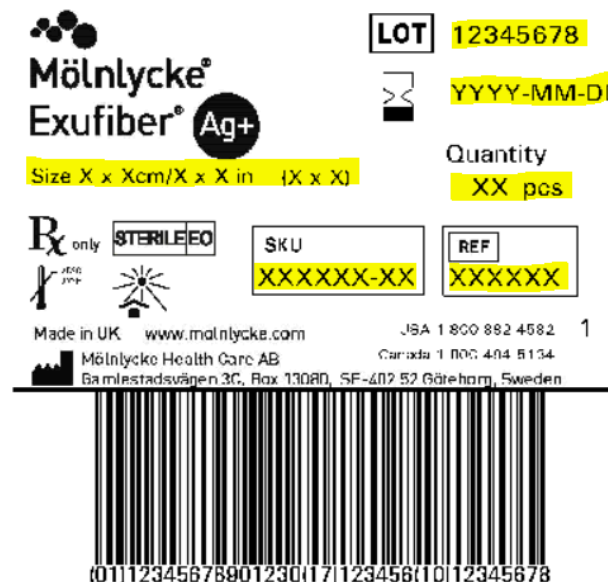

Title: Labelling Specification for Exufiber Ag

Page 7(7)

Example of TRP Label Positioning

