


Infusion - L86 - vented sets

L86 (closed cap, SELF-LOCK)



1. The special design of the spike enables easy penetration of container ports and offers a safe connection to the container. Vented version. Air vent with hydrophobic, sterile air filter and ventilation cap. Standard drip tube (20 drops/ml). The protection cap maintains sterility of the set.
2. The flexible, transparent non-reflective material of the drip chamber with level ring mark enables easy priming and control of drip rate. Equipped with 15 micron fluid filter.
3. Transparent, flexible PVC tubing with plasticiser TOTM 3,0/4,1 mm (ID/OD) in the set (regular bore).
4. Precision roller clamp enables exact adjustment of the desired flow rate maintaining it with minimum variation. Additionally, it is equipped with a clamp/hanging device for tubing and for protection against stick accidents after use of the set.
5. Rotating Male Luer-Lock adapter with SELF-LOCK feature for a safe and comfortable connection. FlowStop protection cap for Male Luer-Lock adapter is equipped with an integrated hydrophobic membrane, ensuring internal sterility prior to use. During priming the FlowStop protection cap remains attached to the Male Luer-Lock adapter and allows air to pass as soon as solution reaches the membrane, the flow stops automatically.

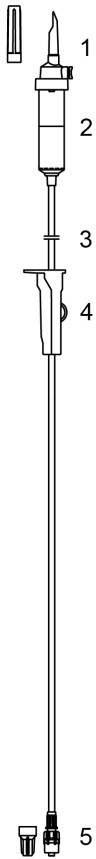

Disclaimer:

The above is indicated only for information. This document is not a specification. The values given are subject to the usual fluctuations. Please follow the instructions for use by the medication manufacturer and the actual SPC (Summary of Product Characteristics). Common hygiene measures and working methods are to be followed at all times. Subject to alterations without prior notice!

REF	Product description	Packaging	GTIN	Filling/Storage volume
43.4535	L86 closed cap, 175 cm PVC (TOTM) tubing 3,0/4,1 mm (ID/OD), SELF-LOCK, FlowStop	1 100 3600	04055698001675 04055698001682 -	7 ml/1 m (40°C) 9 ml/1 m (40°C/2 bar)